


Tandem Repeats Finding

To find tandem repeats, select the *Analyze Find tandems...* context menu item in the *Sequence View* window.

In the opened dialog you can specify the tandem search parameters, the region to search in and the result parameters:

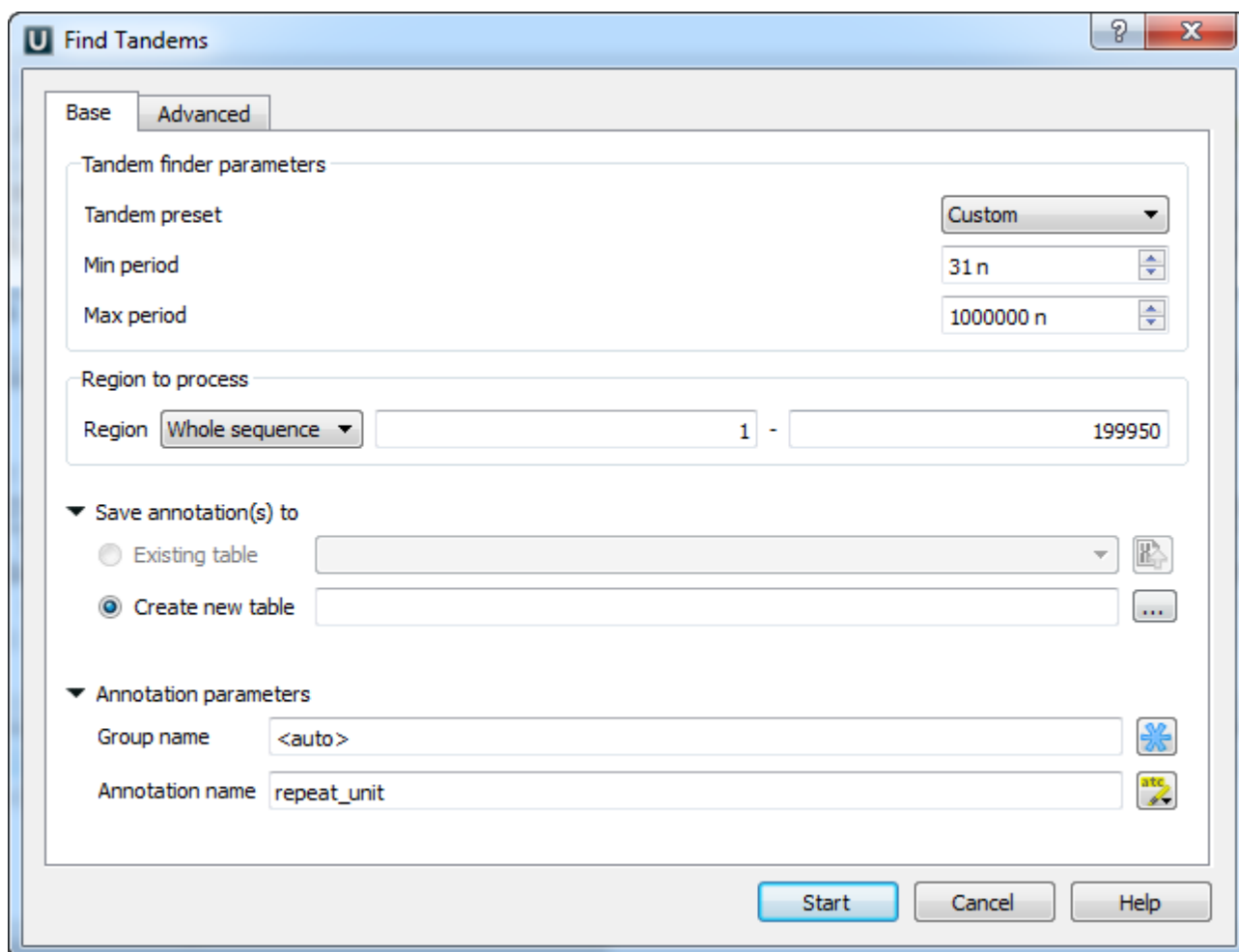


The **Find Tandems** dialog box is shown with the **Base** tab selected. It contains the following sections:

- Tandem finder parameters:**
 - Tandem preset:** A dropdown menu set to **Custom**.
 - Min period:** A text box containing **31 n**.
 - Max period:** A text box containing **1000000 n**.
- Region to process:**
 - Region:** A dropdown menu set to **Whole sequence**.
 - Start/End:** Text boxes containing **1** and **199950**.
- Save annotation(s) to:**
 - ☐ Existing table: A dropdown menu.
 - ☒ Create new table: A text box.
- Annotation parameters:**
 - Group name:** A text box containing **<auto>**.
 - Annotation name:** A text box containing **repeat_unit**.

At the bottom are **Start**, **Cancel**, and **Help** buttons.

The dialog parameters:

Tandem preset — specify the tandem repeats parameters with predefined values by selecting the available preset:

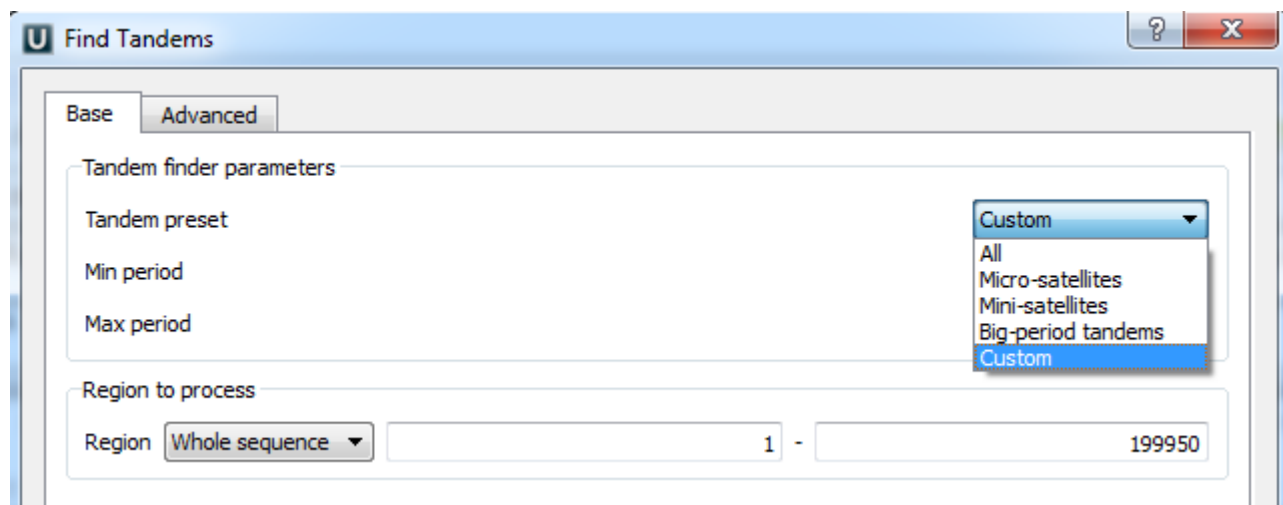

The **Find Tandems** dialog box is shown with the **Base** tab selected. The **Tandem preset** dropdown menu is open, showing the following options:

- Custom
- All
- Micro-satellites
- Mini-satellites
- Big-period tandems
- Custom

The other parameters remain the same as in the first image.

Min period, Max period — the minimum and maximum acceptable repeat length measured in base symbols.

Region to process — specify the region to search in the whole sequence, a custom region or the region of the current selection (if any).

Save annotation(s) to — specify the existing or new annotations table file to save the resulting annotations into.

Annotation parameters — you can change the default group name and annotation(s) name values of the resulting annotation(s).



The screenshot shows a software window with two tabs: 'Base' and 'Advanced'. The 'Advanced' tab is selected. Below the tabs, there is a section titled 'Advanced parameters'. This section contains several settings:

- Repeats identity**: A text box with '100%' and a small up/down arrow icon on the right.
- Algorithm**: A dropdown menu showing 'Suffix array (optimized)' with a downward arrow icon on the right.
- Minimum tandem size**: A text box with '9' and a small up/down arrow icon on the right.
- Minimum repeat count**: A text box with 'x3' and a small up/down arrow icon on the right.
- Show overlapped tandems**: An unchecked checkbox.

Additional search options can be found in the *Advanced* tab:

Algorithm — the algorithm parameter allows to select the search algorithm. The default and a fast one is optimized suffix array algorithm.

Minimum tandem size — the minimum tandem size sets the limit on minimum acceptable length of the tandem, i.e. the minimum total repeats length of the searched tandem.

Minimum repeat count — the minimum number of repeats of a searched tandem.

Show overlapped tandems — check if the plugin should search for the overlapped tandems, otherwise keep unchecked.

- [Tandem Repeats Search Result](#)