


# Circular View Settings

To configure circular view settings go to the *Circular View Settings* tab in the *Options Panel*.

Activate the circular view for a sequence and the following settings will appear:

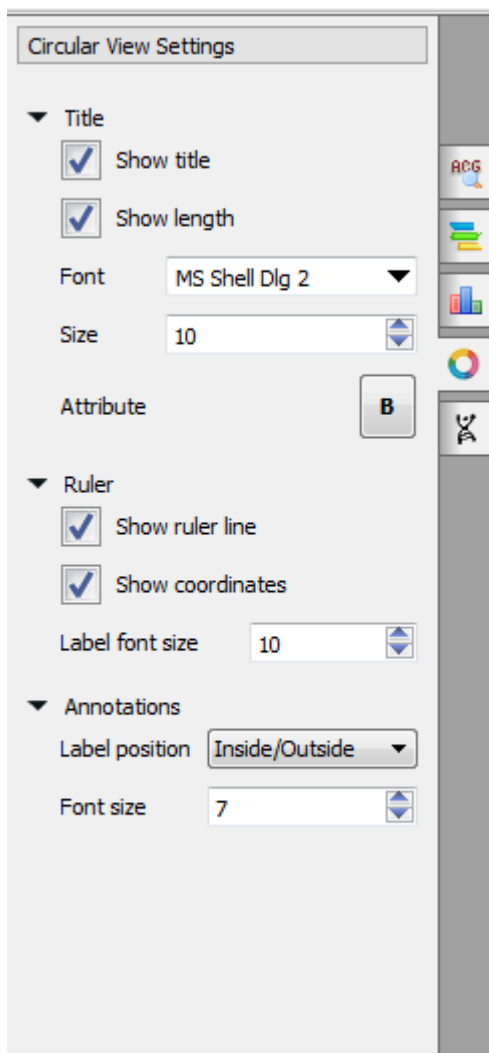

The image shows a software dialog box titled "Circular View Settings". It is organized into three main sections, each with a dropdown arrow to its left:

- Title**: Contains two checked checkboxes for "Show title" and "Show length". Below these are controls for "Font" (set to "MS Shell Dlg 2"), "Size" (set to "10"), and "Attribute" (with a button labeled "B").
- Ruler**: Contains two checked checkboxes for "Show ruler line" and "Show coordinates". Below these is a "Label font size" control set to "10".
- Annotations**: Contains a "Label position" dropdown menu set to "Inside/Outside" and a "Font size" control set to "7".

On the right side of the dialog, there is a vertical toolbar with icons for "ACG", a document, a bar chart, a color wheel, and a molecular structure icon.

In the title section you can show or hide title and length, change font, size and attribute.

In the ruler section you can show or hide ruler line and coordinates and change the label font size.

In the annotation section you can select the label position and change the label size.

The following label positions are available:

- inside - all labels are inside of the annotations
- outside - all labels are outside of the annotations
- inside/outside - if the label can fit the annotation and it is not auto-annotation, it's located inside. Otherwise outside.
- none - no labels at all