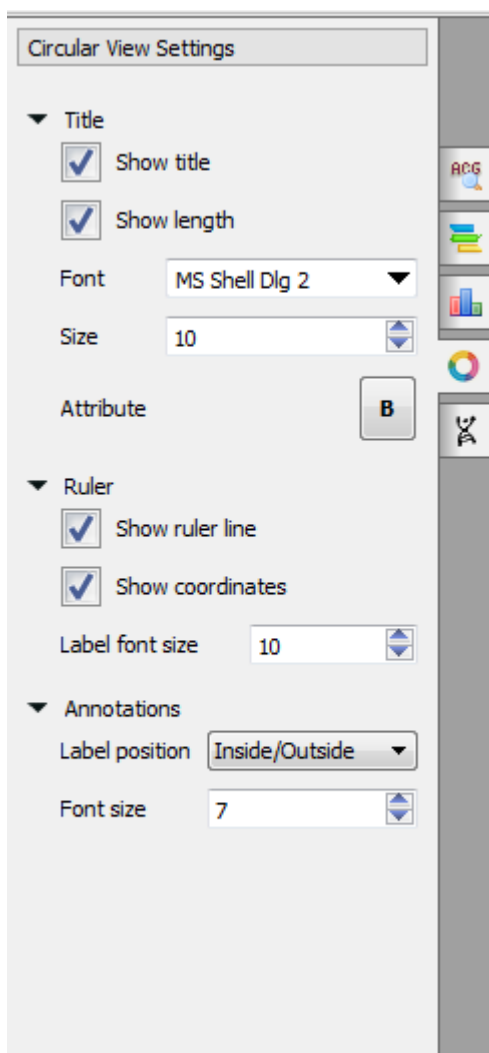


Circular View Settings

To configure circular view settings go to the *Circular View Settings* tab in the *Options Panel*.

Activate the circular view for a sequence and the following settings will appear:

The image shows a software dialog box titled "Circular View Settings". It is organized into three main sections, each with a collapsed arrow icon to its left. The "Title" section contains two checked checkboxes for "Show title" and "Show length", followed by a "Font" dropdown menu set to "MS Shell Dlg 2", a "Size" spinner set to 10, and an "Attribute" button with a bold "B" icon. The "Ruler" section contains two checked checkboxes for "Show ruler line" and "Show coordinates", followed by a "Label font size" spinner set to 10. The "Annotations" section contains a "Label position" dropdown menu set to "Inside/Outside" and a "Font size" spinner set to 7. To the right of the settings area is a vertical toolbar with icons for "ACG" (a magnifying glass over a sequence logo), a sequence logo, a bar chart, a color wheel, and a DNA double helix.

In the title section you can show or hide title and length, change font, size and attribute.

In the ruler section you can show or hide ruler line and coordinates and change the label font size.

In the annotation section you can select the label position and change the label size.

The following label positions are available:

- inside - all labels are inside of the annotations
- outside - all labels are outside of the annotations
- inside/outside - if the label can fit the annotation and it is not auto-annotation, it's located inside. Otherwise outside.
- none - no labels at all